

SORDINA

technology serving life



Sordina Medical Technologies S.r.l.

Tel +39 0444 233750 - Fax +39 0444 233790
info@sordina.com - www.sordina.com

Registered Office

Galleria Pozzo Rosso 13
36100 Vicenza, Italy
REA VI-357804

ST9
Operating Table

Our electro hydraulic
most advances OT
for high – end surgeries



Since more than 140 years Sordina provides its know how to Medical Teams and they're Patients offering high quality and high reliability Operating Tables designed to satisfy numerous needs.

Italian passion, global vision: according to this concept Sordina designed a new generation of high-performances Operating Tables able to fulfill User's expectations. Sordina engineered the new Operating Table line together with Operating Room Teams to ensure product's highest reliability and use friendly.

High reliability, superb movement precision, design, practicality are some of the ingredients Sordina's Team used to develop the new top of the range Operating Table, the ST9 model. The Operating Table ST9 contains the historical and widely confirmed experience adapted to the new Surgical Team needs.

ST9 has an electro-hydraulic heart capable to guarantee high movement reproducibility during Operating Table lifetime; its movements will be always precise and smooth.

Technology serving life, the new generation



IR remote control



Double joint head section



Console equipped with membrane keyboard



Ergonomic base designed to facilitate Surgical Team operations


Technology serving life, the new generation

Antistatic, radiotransparent and latex free cushions are available in memory foam with anti-decubitus effect

To ensure an optimal radiolucency translation movement is available for the whole operating plane

ST9 is a mobile Operating Table with fixed operating plane able to effect a great range of movements optimizing surgical site exposure guaranteeing top class safety and comfort.
The Operating Table it's made by a central module which can be integrated with an ample range of configurations for head, back and leg sections; thanks to this modularity the ST9 it's able to satisfy Patient's positions according to the required surgical specialty.
Sordina is equipping its new ST9 Operating Table line with both wire and IR remote control and a redundant control system installed on board able to execute the whole OT movement range.

Each single section is adjustable independently from the others

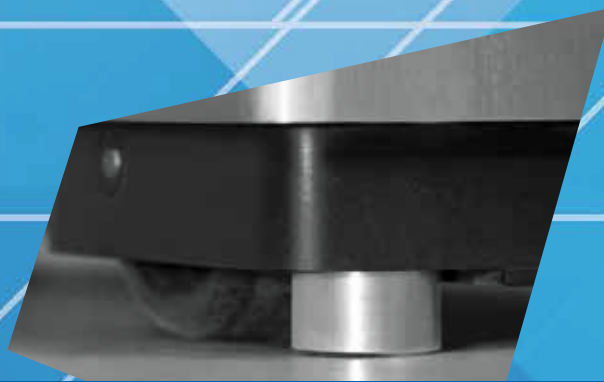

Electro-hydraulic floor locking system




Lock-unlock system



The whole operating plane can be leveled using a single button

Technology serving life, the new generation

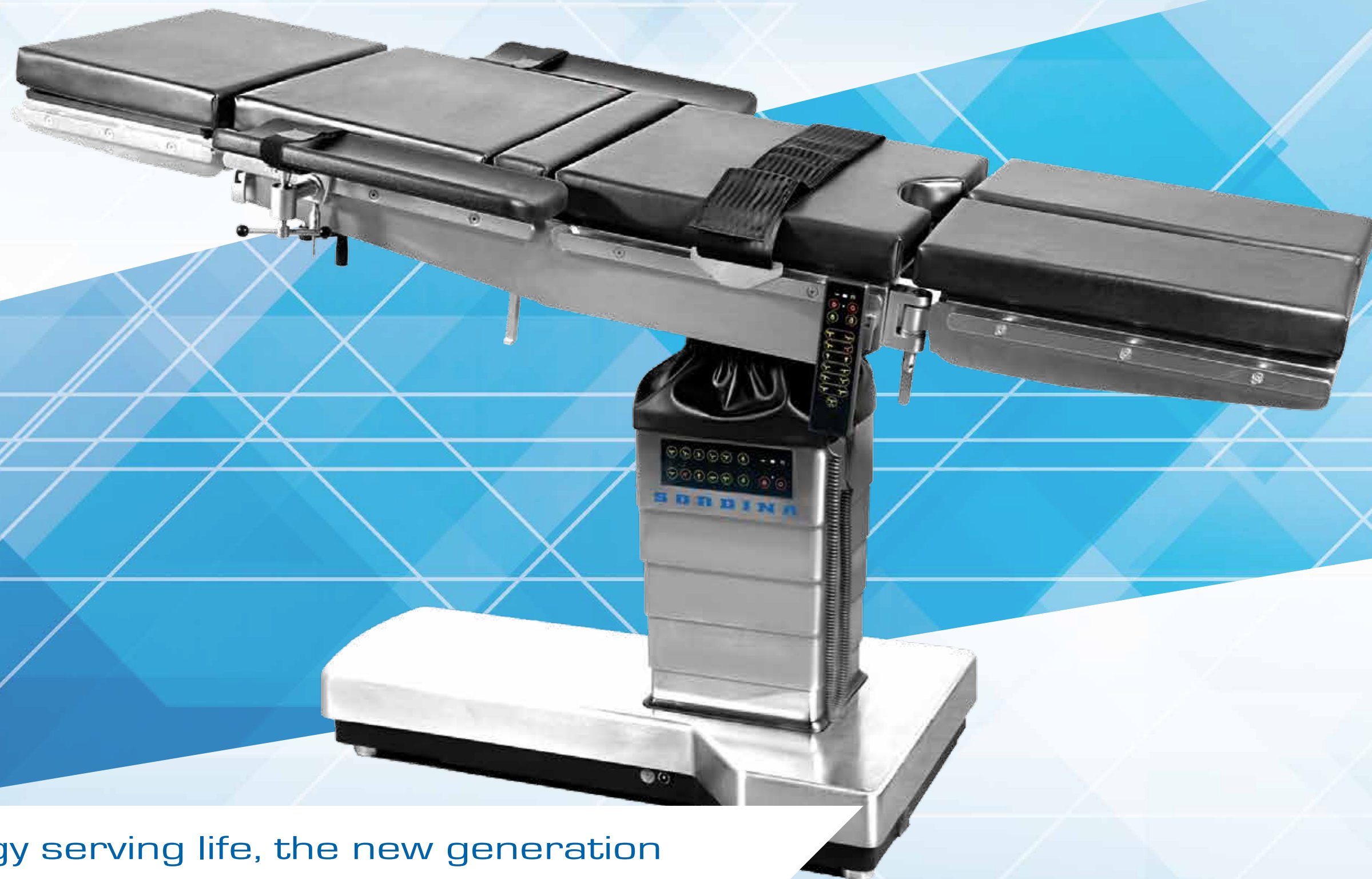

Technology serving life, the new generation

TECHNICAL SPECIFICATION


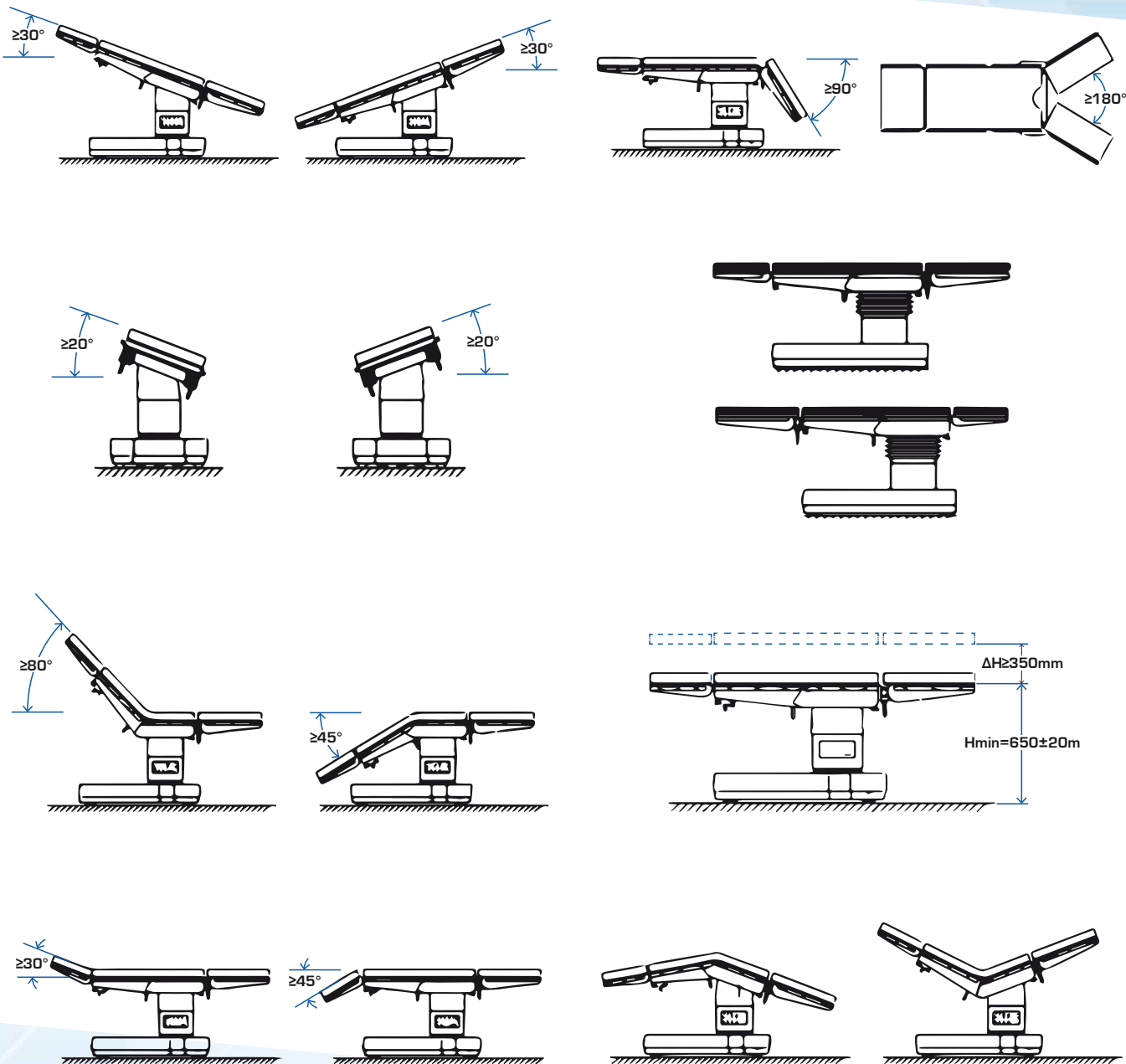




Table Length	2020 mm
Table Width (without side rails)	520 mm
Table Height (without-pads)	650 mm-1000 mm
Trendelenburg/ Reverse	+30°/-30°
Lateral Tilt (left/ right)	20°/20°
Back Plate Up/ Down	+80°/-45°
Head Plate Up/ Down (detachable)	+30°/-45°
Leg Plates Down	90°
Leg Plate Split (detachable)	180°
Kidney-bridge Elevator UP	128 mm
Flex/ Reflex	240°/120°
Loading Capacity	297 Kg
Power Input (Hz)	AC 90 V-264V (50-60)
Cell Output /(Ah)	DC 24V (12)
Longitudinal Sliding	300 mm

Technology serving life, the new generation

SUPERB VERSATILITY

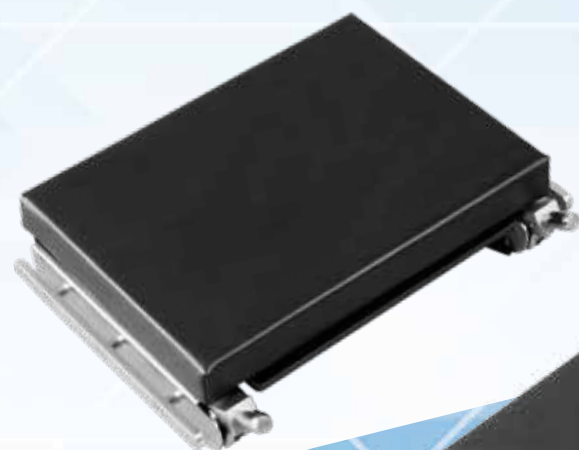

Leg sections allows multiple configurations

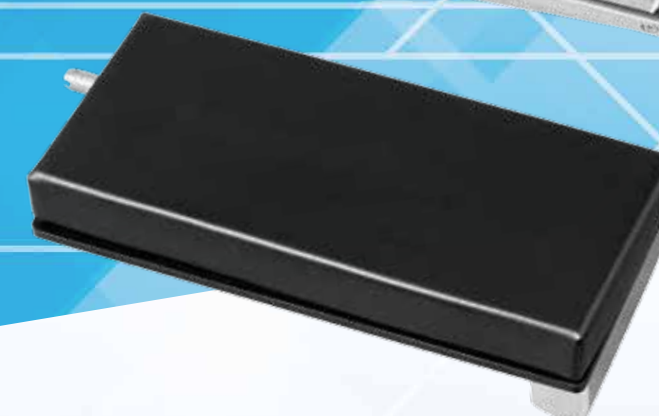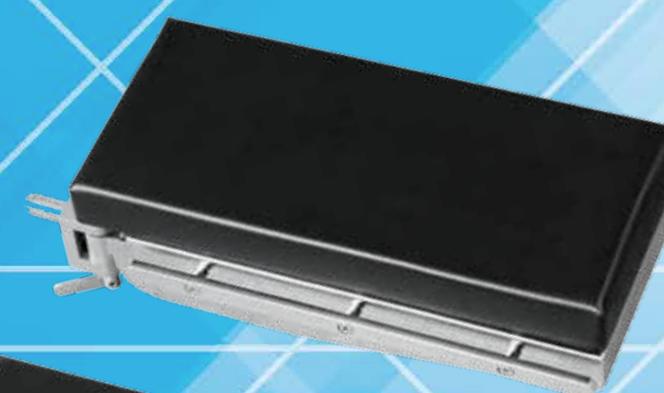

Leg sections are actioned manually by gas spring

Double joint head sections are actioned manually by gas spring

Technology serving life, the new generation



abdominal surgery



spine surgery



thorax surgery



neurosurgery in prone position



gynecological/urological surgery



orthopedic surgery



orthopedic surgery



orthopedic surgery lateral position



knee surgery



neurosurgery

**GUARANTEED AND CERTIFIED
QUALITY SINCE 1875**

SERVICE ASSISTANCE

Numero Assistenza
+39 0444 233 750

service@sordina.com

